

VWR® ERGONOMIC CHAIRS

*When
innovation
becomes
revolution*



5 year Warranty

Ergonomics enters a new era



Lab stool with castors
132-1099



Lab chair with lumbar support
and castors
132-1104



Lab chair with lumbar support,
foot ring and glides
132-1105



Lab stool with foot ring
and glides
132-1100


Ergonomic polyurethane chairs

The seats are made of soft polyurethane with 40 mm thickness and offering exceptional seating comfort.

- Resistant to sparks, sharp tools, chemicals and weak acids
- Easy to clean
- No perspiration effect
- Appropriate for any environment (labs, industries, hospitals)



Scientifically proven benefits


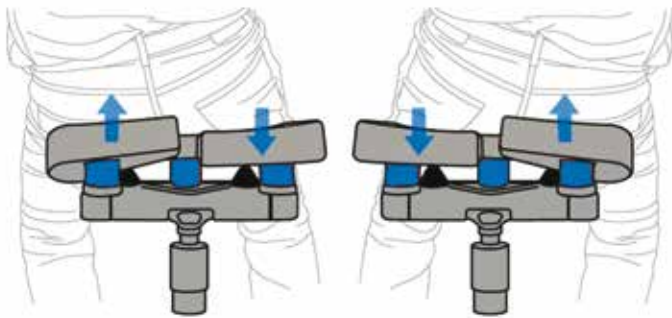


Nobody has a perfectly symmetric morphology. This is why the majority of people have back problems, sometimes with a feeling of heavy legs. The answer to this problem can be found in the main features of this seat.

Above all, a good working position is a dynamic posture, making it possible to move, because moving ensures good bloodstream. The features of the seat help to rebalance the body naturally and improve blood flow in the legs: It consists of a short seat tilted forward, split into two parts, based on silent-blocks forming a triangle.

This seat has a great deformation capacity which provides better distribution of pressure when moving: For instance, when pressing a pedal, tests have shown that the difference in pressure between each leg decreases by 82% compared to a traditional seat, and the pressure under the thigh, which controls the pedal, is reduced by 49%.

This range was developed in close partnership with a company specialising in human motion metrology. Many tests were carried out with different types of seating, people, movements and postures, based on motion capture technology. They helped to highlight the benefits of this seat on posture, and more specifically, on the behaviour of the vertebral rachis.

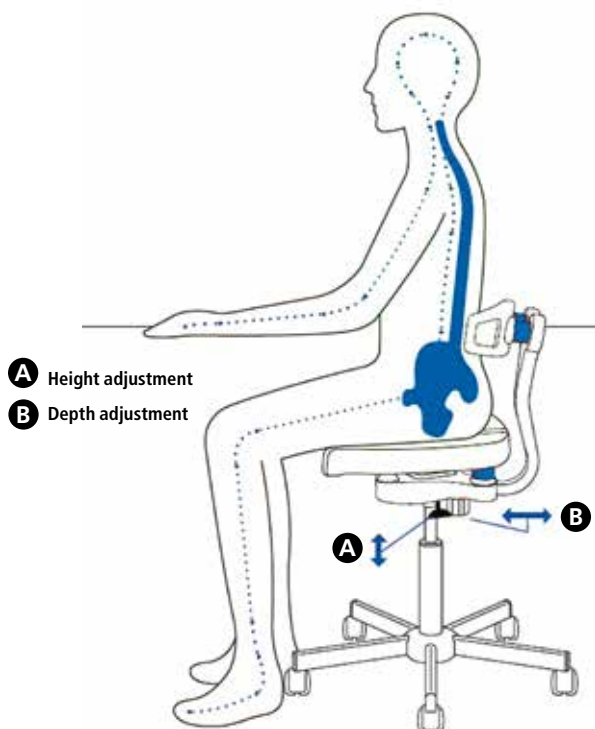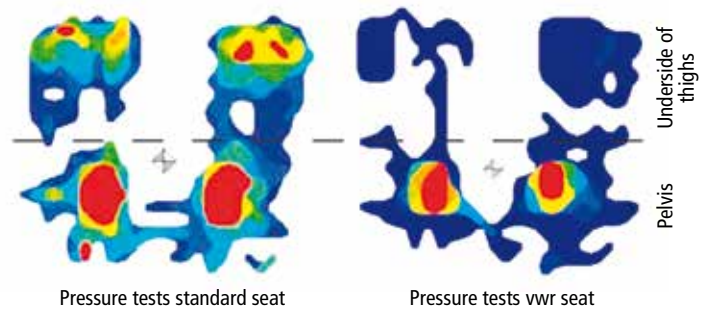

On a work seat, the backrest is only slightly used, or misused. In the latter case, it will become your worst enemy. For this reason, this chair is not equipped with a standard backrest but with a lumbar support, which aims to help you immediately find the perfect posture.

This VWR seat is adjustable in height with a gas-lift column, which is the only manual adjustment, all other adjustments being intuitive, in order to have a seat very easy to use.

The inclination of the seat helps to concentrate supports on the ischiums, which is the anatomical region assigned to this function. The orientation of the seat also facilitates an upright posture and minimises spinal stress. The centre of gravity moves back to focus support under the pelvis and enables a more ergonomic posture.

| Description | Height (mm) | Cat. No. |
|-----------------------------------------------------------|-------------|----------|
| Lab stool with castors for hard floors | 510 - 645 | 132-1099 |
| Lab stool with foot ring and glides | 610 - 855 | 132-1100 |
| Lab chair with lumbar support and castors for hard floors | 510 - 645 | 132-1104 |
| Lab chair with lumbar support, foot ring and glides | 610 - 855 | 132-1105 |

Complies with the main security requirements of:

- NF EN1335-3
- NF D655761
- DIN 68877
- EN 1021-1-2

Austria

VWR International GmbH
Graumannsgasse 7
1150 Vienna
Tel.: +43 1 97 002 0
Fax: +43 1 97 002 600
Email: info@at.vwr.com

Belgium

VWR International bvba
Researchpark Haasrode 2020
Geldenaaksebaan 464
3001 Leuven
Tel.: 016 385 011
Fax: 016 385 385
Email: vwrbe@be.vwr.com

Czech Republic

VWR International s. r. o.
Veetee Business Park
Pražská 442
CZ - 281 67 Stříbrná Skalce
Tel.: +420 321 570 321
Fax: +420 321 570 320
Email: info@cz.vwr.com

Denmark

VWR - Bie & Berntsen
Transformervej 8
2860 Søborg
Tel.: 43 86 87 88
Fax: 43 86 87 90
Email: info@dk.vwr.com

Finland

VWR International Oy
Valimotie 9
00380 Helsinki
Tel.: 09 80 45 51
Fax: 09 80 45 52 00
Email: info@fi.vwr.com

France

VWR International S.A.S.
Le Périgares – Bâtiment B
201, rue Carnot
94126 Fontenay-sous-Bois cedex
Tel.: 0 825 02 30 30 (0,18 € TTC/min)
Fax: 0 825 02 30 35 (0,18 € TTC/min)
Email: info@fr.vwr.com

Germany

VWR International GmbH
Hilpertstraße 20a
D - 64295 Darmstadt
Freecall: 0800 702 00 07
Fax: 0180 570 22 22*
Email: info@de.vwr.com
*0,14 €/Min. aus d. dt. Festnetz

Hungary

VWR International Kft.
Simon László u. 4.
4034 Debrecen
Tel.: (52) 521-130
Fax: (52) 470-069
Email: info@hu.vwr.com

Ireland / Northern Ireland

VWR International Ltd /
VWR International (Northern Ireland) Ltd
Orion Business Campus
Northwest Business Park
Ballycoolin
Dublin 15
Tel.: 01 88 22 222
Fax: 01 88 22 333
Email: sales@ie.vwr.com

Italy

VWR International PBI S.r.l.
Via San Giusto 85
20153 Milano (MI)
Tel.: 02-3320311/02-487791
Fax: 800 152999/02-40090010
Email: info@it.vwr.com

The Netherlands

VWR International B.V.
Postbus 8198
1005 AD Amsterdam
Tel.: 020 4808 400
Fax: 020 4808 480
Email: info@nl.vwr.com

Norway

VWR International AS
Haavard Martinsens vei 30
0978 Oslo
Tel.: 02290
Fax: 815 00 940
Email: info@no.vwr.com

Poland

VWR International Sp. z o.o.
Limkowa 5
80-175 Gdansk
Tel.: 058 32 38 200
Fax: 058 32 38 205
Email: info@pl.vwr.com

Portugal

VWR International -
Material de Laboratório, Lda
Edifício Neopark
Av. Tomás Ribeiro, 43- 3 D
2790-221 Carnaxide
Tel.: 21 3600 770
Fax: 21 3600 798/9
Email: info@pt.vwr.com

Spain

VWR International Eurolab S.L.
C/ Tecnología 5-17
A-7 Llinars Park
08450 - Llinars del Vallès
Barcelona
Tel.: 902 222 897
Fax: 902 430 657
Email: info@es.vwr.com

Sweden

VWR International AB
Fagerstagatan 18a
163 94 Stockholm
Tel.: 08 621 34 00
Fax: 08 621 34 66
Email: kundservice@se.vwr.com

Switzerland

VWR International GmbH
Lerzenstrasse 16/18
8953 Dietikon
Tel.: 044 745 13 13
Fax: 044 745 13 10
Email: info@ch.vwr.com

Turkey

VWR International Laboratuvar
Teknolojileri Ltd.Şti.
Orta Mah. Cemal Gürsel Caddesi
34896 Pendik - Istanbul
Tel.: +90216 598 2900
Fax: +90216 598 2907
Email: info@tr.vwr.com

UK

VWR International Ltd
Customer Service Centre
Hunter Boulevard - Magna Park
Lutterworth
Leicestershire
LE17 4XN
Tel.: 0800 22 33 44
Fax: 01455 55 85 86
Email: uksales@uk.vwr.com

Australia

VWR International, Pty Ltd.
Level 1, Unit 1a/60 Enterprise Place
Tingalpa, Queensland, 4173
Tel.: 1300 727 696
Fax: 1300 135 123
Email: sales@au.vwr.com

China

VWR International China Co., Ltd.
Shanghai Branch
Room 256, No. 3058 Pusan Road
Pudong New District
Shanghai 200123
Tel.: +86-21-5898 6888
Fax: +86-21-5855 8801
Email: info_china@vwr.com

India

VWR Lab Products Private Limited
No.139, BDA Industrial Suburb,
6th Main, Tumkur Road, Peenya Post,
Bangalore, India – 560058
Tel.: +91-80-28078400
Fax: +91-80-28078410
Email: vwr_india@vwr.com

New Zealand

Global Science - A VWR Company
241 Bush Road
Albany 0632, Auckland
Tel.: 0800 734 100
Fax: 0800 999 002
Email: sales@globalscience.co.nz

Singapore

VWR Singapore Pte Ltd
18 Gul Drive
Singapore 629468
Tel: +65 6505 0760
Fax: +65 6264 3780
Email: sales@sg.vwr.com

GO TO **VWR.COM** FOR
THE LATEST NEWS,
SPECIAL OFFERS AND
DETAILS OF YOUR LOCAL
VWR DISTRIBUTOR